


## Advanced to a Higher Quality and Data Rate!!

# Smart LDV III Model 8743/8743-S

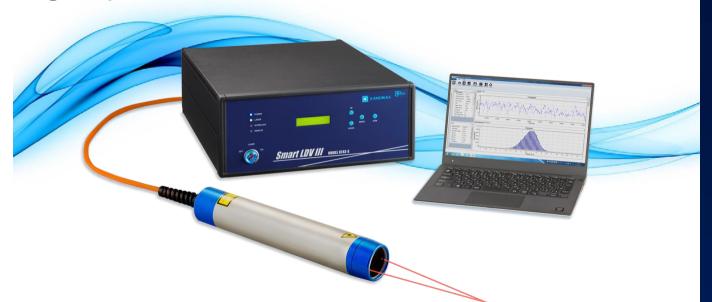


## Features:

- Enhanced receiving sensitivity enabling high data-rate measurement
- Up to 60,000 velocity data/sec.
- High-speed data transfer by USB3.0
- Probe designed as all-in-one, no alignment necessary

## **Applications:**

- Aerodynamic and hydraulic property measurement
- Measurement requiring high time-resolution
- Comparison with CFD
- PIV accuracy tests

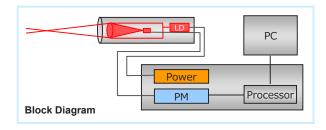

## **Specifications**

•	
Flow velocity range	-40m/s~260m/s
Ontical System	(f=400mm, Model 8743-S)
Optical System  Laser	LD λ=660nm, 60mW
Lasei	150mm, 200mm, 250mm, 300mm,
Focal length	350mm, 400mm
Measurement volume size	0.13mm x 1.3mm (f=200mm)
Measurement method	Back scatter / Forward scatter(Option)
Probe size	Dia. 61mm x 345mm
Shift frequency	Model 8743: Without Frequency Shifter Model 8743-S: 0.01 – 10MHz
Power supply	AC100 - 240V
Signal Processor	
Signal processing	8bit FFT (512,256,128point)
Frequency band	1kHz~40MHz (8 ranges)
Max data rate	60,000 speed data/sec*
Validation	Burst spectrum ratio
Interface	USB3.0
Software	
Max. number of data	100,000
Real time monitor	Burst waveform
	Burst spectrum
	Velocity histogram
Analysis function	Mean flow velocity, Turbulent intensity,
	Skewness factor, Flatness factor
	Velocity histogram, Time-series display
Data output	CSV format
Supported Operating System	Windows 7 / 8 / 10 (64bit only)
	Japanese / English
	*Depending on measurement condition

KANOMAX USA, INC. 219 US HWY 206 ANDOVER, NJ 07821 USA **TEL:** 973-786-6386 FAX: 973-786-7586

E-MAIL: SALES@KANOMAX-USA.COM **URL: WWW.KANOMAX-USA.COM** 

#### **DISTRIBUTED BY:**



## **Options**

## **Traverse System**

**Automated Traverse System for positioning the optical system** 

Easy measurements without the hassle of changing the measurement location manually.

- Automatic measurements from the LDV software
- Moves between each measurement point with high positioning accuracy
- Operable also in manual mode

 $\begin{array}{ccc} \text{Moving axis} & & X,\,Y,\,Z \text{ axis} & ^* \\ \text{Stroke} & & 500 \text{ mm} & ^* \end{array}$ 

Positioning accuracy ±0.02 mm (X axis) \*\*
Drive system Stepping motor
Controller LDV Software

(traverse-compatible version)

\* Reference examples



## **Corner Cube Mirror**

## Corner cube Mirror for better SNR of data


Designed to be placed at the opposite side of the laser to reflect strong front scattering light to the optical receiver in order to improve the signal-to-noise ratio.

Focal length 200 mm, 250 mm, 300 mm, 350 mm, 400 mm

Effective diameter Dia. 50 mm

## **Application Example: Cylinder Wake Measurement**

## **Instrument Set-up**

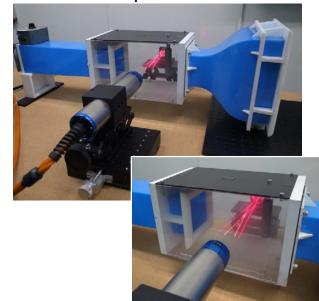


#### **Measurement Condition**

Tracer particle: Approx. dia. 5.0 µm

Focal length of probe: 200 mm
Frequency shifter: Available
Measurement point: Cylinder wake

#### Result



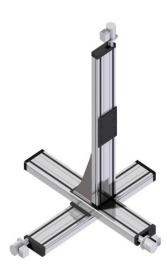

<sup>\*\*</sup> Positioning accuracy for Y and Z axis varies depending on the load